




The **NIKE** couch is a professional couch used by rehabilitation and physiotherapy professionals for examination and Postural Re-Education treatments.



Frame color RAL 9001

TECHNICAL SPECIFICATIONS

| | |
|---|----------------------|
| Sections | 5 |
| Lateral offset | 17 cm |
| Width | 65 cm |
| Length | 190 cm |
| Max working load, evenly distributed | 200 kg |
| Height – electrically or hydraulically adjustable | 45 - 90 cm |
| Middle section angle | with gas spring +20° |
| Footrest angle | with gas spring +75° |
| Trendelenburg position | + 20° |
| Power supply | 230 V 50/60 Hz 1,5A |
| Actuator thrust | 8000 N |
| IP Code | IP 66 |
| Headrest angle – with gas spring | +60° – 20° |
| Breathing hole and bung | as standard |
| Removable armrests rotating by | 180° |
| Weight | 70 kg |

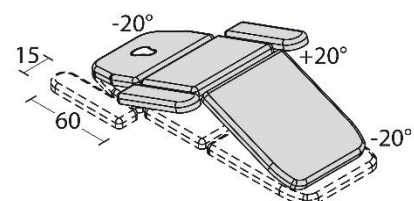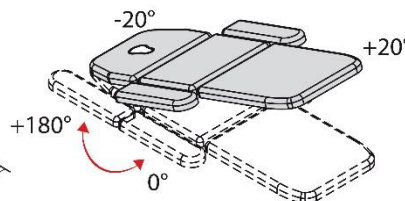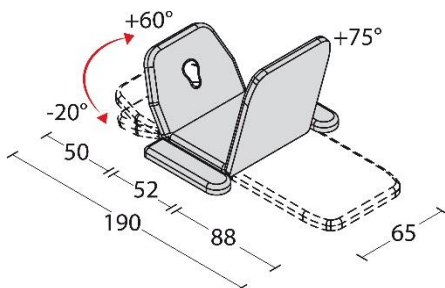

FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec :

fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

NIKE COUCH MODELS

| Description | Code |
|--|------------|
| Nike with foot pedal. | 1117C10P1S |
| Nike with all around foot bar switch . | 1117C10C1S |
| Hydraulic Nike. | 1117C10R0S |

HEIGHT ADJUSTMENT SYSTEMS

NIKE with 24V foot pedal



NIKE with all around foot bar switch



NIKE – hydraulic



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec :
fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

TECHNICAL SPECIFICATIONS – PADDING

| | |
|------------------------|----------------------|
| Name of material | Polyurethane foam |
| Reaction to fire class | 1 IM (highest) |
| Density | 30 kg/m ³ |
| Height | 40 mm |



AVAILABLE UPHOLSTERY COLOURS



Beige P1



Light grey G2



Taupe P2



Dark grey G3



Yellow G1



Mustard G4



Orange A1



Red R1



Burgundy B2



Violet V4



Dark blue B3




Blue B1



Dark brown M1



Earth M2

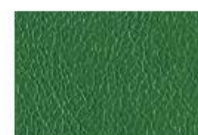

Green V1



Forest green V2



Leaf green V3



Black N1

FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec :

fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

TECHNICAL SPECIFICATIONS – UPHOLSTERY

| FEATURE | TEST METHOD | UNIT | LIMIT VALUE |
|--|--|---------------|--|
| Tensile strength | UNI EN ISO 1421 Method 1 | N | Lengthwise ≥ 210 Crosswise ≥ 90 |
| Elongation strength | UNI EN ISO 1421 Method 1 | % | Lengthwise ≥ 85 Crosswise ≥ 140 |
| Determination of tear force | EN ISO 13937-2: 2000 | daN | Lengthwise $\geq 2,5$ Crosswise $\geq 2,2$ |
| Colour fastness to crocking | UNI EN ISO 105 X12: 2003 UNI EN ISO 105 D02: 1998 | grey scale | Dry 4 Bleach 4/5 Ethyl alcohol 3/4 |
| Determination of abrasion resistance: Martindale abrader | UNI EN ISO 5470-2: 2004 | cycles | Method 1 (direct) 65.000 Method 2 (reverse) 210.000 |
| Reaction ignition Fireproof class | UNI 9175: 1987 UNI 9175/FA1: 1994 | | CLASSE 1.IM |
| Free and hydrolyzed formaldehyde | UNI EN ISO 14184-1: 2011 | mg/Kg | < 3 |
| Phtalates | UNI EN ISO 14389: 2014 | % | < 0,06 |
| Arylamines | UNI EN ISO 14362-1: 2012 | mg/Kg | Not detectable |
| Cleaning | Clean with damp soapy cloth and rinse well with clean water. Do not use solvents, bleaches, synthetic detergents, polishes or aerosol spray. | | |
| Specific features | Antibacterial Fingicide | | |
| Composition | Plasticised PVC 88% - Polyester 12% | | |

FIRE RESISTANCE ACCORDING TO:
IMO - 2010 FTP CODE Part 8
"Test for upholstered furniture"
par 7.2 SMOULDERING CIGARETTE TEST
par 7.3 PROPANE FLAME TEST
EN 1021-1: 2006 smouldering cigarette test
EN 1021-2: 2006 match flame equivalent test

FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec :

fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

NIKE COUCH POSTURAL SERIES OPTIONAL ACCESSORIES

Code 1210010 PAPER ROLL HOLDER

Metal structure to be placed under the final section of the couch.



Code 1211010 CASTORS

Four retractable castors to be used only to move the couch.
Ø 60 mm.



Code 1217010 SAFETY MAGNET

Removing the magnet all electrical regulation functions are prevented. This avoids any accidental activation of the actuators.

Code 1217020 SAFETY DOUBLE CLICK WITH MAGNET

A software that allows the movement of the actuators only through a double input of the control. This avoids any accidental activation of the actuators.



Code 1217030 SWISS PLUG



Code 1217031 UK 3-PIN PLUG



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec :

fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956